

WOMEN'S IMAGING Ultrasound Table

- Full featured 3 section table top design for all Women's Imaging ultrasound modalities
- 550 lb (250kg) bariatric rated lift & load capacity
- Twin tower design for added strength and rigidity
- Sleek base design promotes unobstructed compatibility with ultrasound equipment ensuring better ergonomic conditions for the sonographer



shown in Forest fabric

Inspired by the leading ergonomic experts in sonography



SOUND | ERGONOMICS

- Decrease sonographer injury
- Features minimize exam times



FM584901

made in the USA
with US & imported parts



OAKWORKS®
MEDICAL
Position for Success

WOMEN'S IMAGING Ultrasound Table

Recommended for:

- OB/GYN exams
- All abdominal & superficial structures exams
- Musculoskeletal exams
- Vascular exams
- Biopsies



shown in Opal fabric with optional second Adjustable Back support, Central Base Locking System and Side Rails

WOMEN'S IMAGING TABLE

#66127-T* 27" W #61848-T* 30" W

FEATURES	
3 Section Top	27" (69cm) or 30" (76cm) wide x 72" (183cm) long
Backrest	34" (86cm) long
Seat	24" (61cm) long
Leg Rest	13" (33cm) long
Motion 1	Electric Height Range 22" - 38" (56-97cm)
Motion 2	Electric Trendelenburg/Reverse Tilt $\pm 15^\circ$
Motion 3	Electric Fowler Backrest 0-75°
Motion 4	Manual 4 position Leg Rest 0-90°
Controls	Multi-function hand control
Weight Capacity	550 lb (250kg)
Product Weight	359 lb (163kg)
Padding	3" (8cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	4" (10cm) locking casters, Adjustable access back support (right) 0-85°, Paper roll holder, Stirrups (Retractable foot rests)
Warranty	3 years parts, 2 years labor

OPTIONS must be ordered with table.		#
Accessory Mounting Plate	Fowler end	62123
Adjustable Access/Back Support	Left side 0-85°	73133
Foot Control	Height, Trendelenburg, Fowler	61901
Foot Control	Height, Trendelenburg, Fowler, Cardio Panel	65679
Central Base Locking System	Not available with Patient Lift Compatible Base	61777
Non-USA electric system	please contact us for order details	73134
Patient Lift Compatible Base	Not available with Central Base Locking System	65893
Procedure Tray 30" (76 cm) Width Only	12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters)	73127
Side Rails (set of 2)	23.5" (60cm) long x 6" (15cm) tall (from surface of table)	61781
Sliding Cardio Panel - Electric	30" (76cm) width table, left side	74080
Sliding Cardio Panel - Manual	30" (76cm) width table, left side	73130
ACCESSORIES order any time.		#
Detachable Arm Board	Use with T-rail adapter 5" (13cm) wide x 20.5" (52cm) long	62192
QuickLock™ Head Rest (black) & AeroCel Pad	Comfortable, adjustable prone or supine support (requires mounting plate)	PKG7554*
Table Extender	24" (61cm) wide x 12" (30cm) long (requires mounting plate)	63487

*Fabric color code must be specified.

ULTRASOUND FEATURES & OPTIONS

CENTRAL BASE

LOCKING SYSTEM #61777

Easily lock your table into place with our unique system that allows locking from both sides of the table.



OPTIONAL

SIDE RAILS #61781

Increase stability and security for your patient. Fold down out of the way when not in use. 23.5" (60cm) long x 6" (15cm) tall (from surface of table)



OPTIONAL

STIRRUPS #71001

Adjust to a wide range of angles & lengths to accommodate any patient or procedure. Retractable design stores inside the table.



STANDARD

FULL FUNCTION FOOT CONTROLS

#enquire for part number

You can add a foot control that controls all functions.



Multi-function foot control

PATIENT LIFT COMPATIBLE BASE #65893

This option allows the use of patient lifts and easy wheel chair transfers to minimize or eliminate caregiver injuries associated with patient transfers. Must be ordered with table.

- Includes larger, 5" (13cm) diameter Heavy Duty locking casters
- Provides 6.5" (17cm) clearance for patient lift
- Central Base Locking System not available with this option



OPTIONAL

SLIDING CARDIO PANEL patent pending

Electric: #74080 Manual: #73130

The sliding cardio panel allows you to adjust the opening to the size you need, up to 15" (38cm). 10" (25cm) wide x 28" (71cm) long. Great for left or right handed scanning. Can also be used as an armrest. Available in manual or electric adjustment options. Available on left side of 30" (76cm) width tables.



OPTIONAL

ADJUSTABLE ACCESS BACK SUPPORT

Left: #73133
Right: #73132

Adjustable with 5 positions 0-85°, allowing easy access for sonographer to the far side and back support for patients. 10" (25cm) wide x 15.5" (39cm) long



(RIGHT)

(LEFT)

STANDARD

OPTIONAL

ACCESSORY MOUNTING PLATE #62123

Add a face rest or table extender to your table.



OPTIONAL

PROCEDURE TRAY #73127

Retractable, removable stainless steel. 12.25" (6cm) x 10.25" (26cm) and holds 2 quarts (1.9 liters). Available on 30" (76cm) wide tables.



OPTIONAL

ULTRASOUND ACCESSORIES

DETACHABLE ARM BOARD #62192-T*

Attaches to the T-Rail Adaptor and rotates over 140 degrees. 5" (13cm) wide x 20.5" (52cm) long (Requires T-Rail Adaptor)



T-RAIL ADAPTOR #61902

Multiple locations allow you to choose the best spot for this easily installed & removed accessory capable of handling up to 150 lb (68kg). Can be used with any American standard accessory. 8" (20cm) long

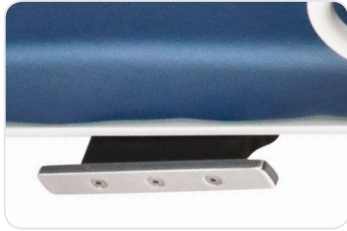


TABLE EXTENDER #63487-T*

Lengthens the table by 12" (31cm) to 84" (213cm) allowing for access to the patient's neck for carotid artery exams and thyroid exams. 24" (61cm) W x 12" (31cm) L



QUICKLOCK™ HEAD REST & AERO-CEL PAD #PKG7554*


The platform allows precise cervical flexion and height positioning over a wide range for both prone and supine treatments. (Requires Accessory Mounting Plate)



TERRATOUGH™ FABRIC

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB117 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.

PADDING

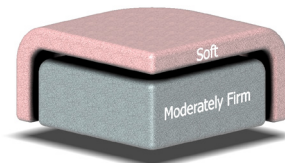

3" Comfort Foam™
More comfort and support



We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

CONTACT:

View our full line of
Medical products

<https://www.oakworksmed.com/pdfdocs/oakworks-medical-catalog/>

717.235.6807 www.oakworksmed.com email: customerservice@oakworks.com