


# PT400

## *Physical Therapy Table*


- 550 lb lifting capacity
- Powered center break
- Can be used as a chair or table
- Top can be put in Trendelenburg position



Shown in Ocean fabric with  
Optional Locking Twin Casters

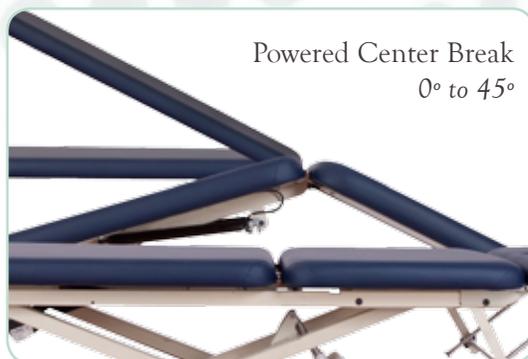

Powered Center Break  
0° to 45°

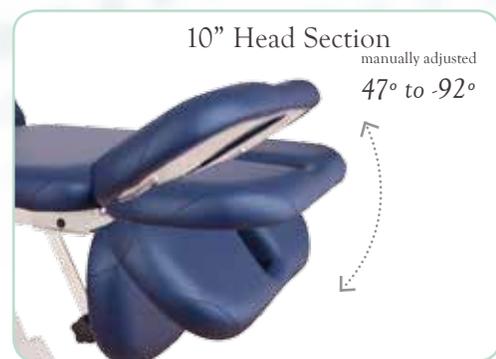


10" Head Section  
manually adjusted  
47° to -92°



ISO 13485:2003  
FM 584901

made in the USA  
with US & imported parts

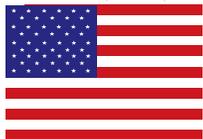

OAKWORKS<sup>®</sup>  
**MEDICAL**  
*Position for Success*

# PT400

#70508-T\* 29" W 1.75" Firm Response    #73393-T\* 29" W 2.5" Comfort  
 #71775-T\* 31" W 1.75" Firm Response    #73392-T\* 31" W 2.5" Comfort

FEATURES		
3 Section Top	29" (74cm) or 31" (79cm) wide x 73" (185cm) long	
Backrest	39" (99cm) long	
Seat	23" (58cm) long	
Head	11" (28cm) long	
Motion 1	Electric Height Range: 17" - 35" (43-89cm)	
Motion 2	Electric Center Break 0° - 45°	
Motion 3	Manual Lift-Assist Backrest 0° - 85°	
Motion 4	Manual 4 position head section 47° to -92°	
Controls	Multi-function foot control	
Weight Capacity	550 lb (250kg)	
Product Weight	198 lb (90 kg)	
Padding	1.75" (4cm) Firm Response or 2.5" (6cm) Comfort Foam™	
Fabric	Premium, PVC free TerraTouch™	
Warranty	3 years parts, 2 years labor	
Safety Listings	FDA listed, CE marked, cETLus approved, CB Certificate	
OPTIONS <i>must be ordered with table.</i>		#
Locking Twin Casters (one end of the table)	To easily move your table in the room	63400
Non-USA electric system	call for ordering details	73145
Second Foot Control	Controls all powered movements	61900
Hand Control	Controls all powered movements	69405

## STIFF TWIN BEAM CONSTRUCTION



Stiffer than cantilever designs and with fewer moving parts, our design supports all your work with a 17" - 35" height range and a powerful ETL listed motor.

## FIRM RESPONSE FOAM



Firm Response foam is designed to deliver a firm yet comfortable surface that maximizes Physical Therapists' ergonomic efficiency during examination and intervention procedures. Additionally, this surface minimizes Manipulation Dry Hits with its super high density, twin layer design. The Firm Response foam enables increased reliability of patient examinations.



### Watch the video!

Dr. Robert DuVall,  
PT, DHS, ATC, OCS, SCS, FAAOMPT CSCS



## TERRATOUGH™ FABRIC

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables and passes testing for Cytotoxicity, Skin Irritation and Sensitization as well as CA TB117 Fire Retardance testing. For disinfecting purposes you may use Protex, MadaCide, MadaWipe, Accell TB, Virox® (Canada), PDI Duper Sani-Cloth Germicidal wipe or a diluted 1/10 bleach solution.

T26 Pure White	T09 Butter	T21 Saffron	T12 Opal	T18 Clay	T28 Taupe	T07 Earth
T27 Espresso	T08 Ruby	T19 Orchid	T05 Heron	T20 Ocean	T17 Sapphire	T04 Sky Blue
T30 Seafoam	T22 Sage	T24 Blue Grass	T14 Forest	T13 Stone	T29 Pewter	T01 Coal

*We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.*

### CONTACT:

### View our full line of Medical products



Scan this code with your smart phone

#72703



717.235.6807    www.oakworksmed.com    email: customerservice@oakworks.com