

# EC SERIES

## *Procedure Chair*

- Powered height, lift-assist backrest & manual leg rest
- 19" to 37" height range
- 550 lbs. weight capacity



Shown in Opal Fabric with  
Optional Arm Rests



FM584901

made in the USA  
with US & imported parts



**OAKWORKS**<sup>®</sup>  
**MEDICAL PRODUCTS**

# EC SERIES *Procedure Chair*

The EC Series Procedure Chair is a cost effective solution for environments where flexibility and lower cost are important. This versatile chair can also be put in a flat table position. A wide range of options and accessories allow you to customize the product to your needs.



Shown with Optional Stirrups, and Side Rails.

The EC Series Chair is a great solution for OB/GYN when configured with stirrups, offering the comfort and flexibility to handle a busy practice.

## EC SERIES PROCEDURE CHAIR #74196-T\* 2.5" Comfort #74380-T\* 4" Comfort

FEATURES	
3 Section Top	29" (74cm) wide x 73" (185cm) long
Back rest	29" (74cm) wide x 39" (99cm) long
Seat	29" (74cm) wide x 23" (58cm) long
Leg rest	29" (74cm) wide x 10" (25cm) long
Leg rest (w/stirrups option)	19" (48cm) wide x 10" (25cm) long
Motion 1	Electric Height Range: 19" - 37" (48-94cm)
Motion 2	Lift-Assist Backrest: 0-90°
Motion 3	Manual Leg Rest: 0-90°
Foot Control	Controls height
Weight Capacity	550 lb (250kg)
Product Weight	165 lb (75kg)
Padding	2.5" (6cm) or 4" (10cm) Comfort Foam™
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, CE marked, cETLus listed


OPTIONS <i>must be ordered with table.</i>		#
Stirrups (retractable foot rests)	Pull out from under seat section, fully adjustable	74201
Side Rails	Drop away (Can not be used with Arm Rests)	63265
Arm Rest	Fold down design (Can not be used with Side Rails)	67648
Locking Twin Casters (head end of table)	To easily move table in & around the room	63400
2nd Foot Control	Height	69615
Hand Control	Controls all powered movements	69405
Non-USA electric system	Please contact us for order details	73134
ACCESSORIES <i>order any time.</i>		#
T-Rail Adapter	Attaches to the backrest for use with USA format T-rail attachable accessories	61902
Quick Lock Head Rest & AeroCel Pad	Comfortable, adjustable prone or supine support	PKG 6914-T*
Paper Roll Holder	Vanilla - Accommodates 21" (53cm) roll (roll not included) - Not FDA listed	11580

\* Fabric color code must be specified.

## STIRRUPS #74201

Adjust to a wide range of angles & lengths to accommodate any patient or procedure. Retractable design stores inside the table. The leg section can be used to support patient's legs while transferring them to the stirrups, then folded down out of the way.

OPTIONAL

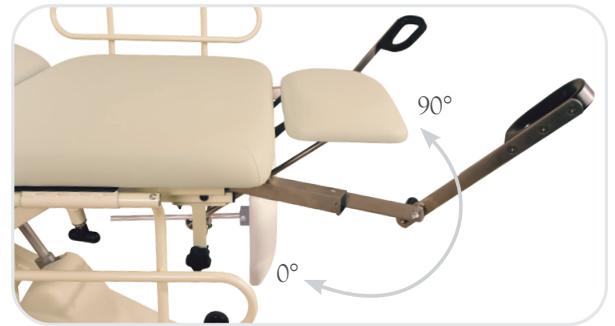

## SIDE RAILS #63265

Increase stability and security for your patient. Fold down out of the way when not in use. 23.5" (60cm) long x 6" (15cm) tall (from surface of table)



OPTIONAL

## ARM REST #67648

Fold down arm rests help the patient feel more secure. They fold down completely out of the way to ensure great access.



OPTIONAL

## LOCKING CASTERS #63400

Easily move your table in and around the room.



OPTIONAL

## T-RAIL ADAPTER #61902

Multiple locations allow you to choose the best spot for this easily installed & removed accessory capable of handling up to 150 lb (68kg). Can be used with any American standard accessory. 8" (20cm) long



ACCESSORY

## QUICKLOCK HEAD REST AND AEROCEL PAD PKG 6914-T\*

This removable head rest is adjustable in angle, height, and distance from the table.



ACCESSORY

## PAPER ROLL HOLDER #11580

This simple design is adjustable to accommodate different size rolls up to 21" (53cm) wide. (not FDA listed)



ACCESSORY

### PADDING



**2.5" Comfort Foam™**  
*Comfort and support*



**4" Comfort Foam™**  
*Superior comfort and support*



Large height range facilitates patient transfer & practitioner ergonomics.

### TERRATOUGH™ FABRIC

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB117 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.



Pure White-T26



Opal-T12



Clay-T18



Earth-T07



Ocean-T20



Sage-T22



Ruby-T08



Orchid-T19



Heron-T05




Coal-T01



Sapphire-T17



Blue Grass-T24

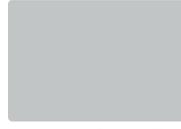

Stone-T13



Pewter-T29

*We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.*

717.759.3146 [www.oakworksmed.com](http://www.oakworksmed.com) email: [customerservice@oakworks.com](mailto:customerservice@oakworks.com)